ABclonal® www.abclonal.com

α-Smooth Muscle Actin (ACTA2) Rabbit mAb

Catalog No.: A17910 Recombinant 129 Publications

Basic Information

Observed MW

42kDa

Calculated MW

42kDa

Category

Primary antibody

Applications

ELISA,WB,IF/ICC

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC0298

Background

This gene encodes one of six different actin proteins. Actins are highly conserved proteins that are involved in cell motility, structure, integrity, and intercellular signaling. The encoded protein is a smooth muscle actin that is involved in vascular contractility and blood pressure homeostasis. Mutations in this gene cause a variety of vascular diseases, such as thoracic aortic disease, coronary artery disease, stroke, and Moyamoya disease, as well as multisystemic smooth muscle dysfunction syndrome.

Recommended Dilutions

WB 1:500 - 1:2000

IF/ICC 1:50 - 1:200

Immunogen Information

Gene ID59

Swiss Prot
P62736

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 1-100 of human α -Smooth Muscle Actin (ACTA2) (P62736).

Synonyms

ACTSA; α-Smooth Muscle Actin (ACTA2)

Contact

<u>a</u>	400-999-6126
\bowtie	cn.market@abclonal.com.cn
$\overline{\Box}$	www.ahclonal.com.cn

Product Information


SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

Validation Data

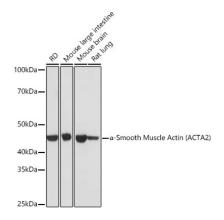

Western blot analysis of various lysates, using α -Smooth Muscle Actin (ACTA2) Rabbit mAb (A17910) at 1:1000 dilution.


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 1s.

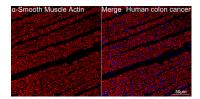

Confocal imaging of paraffin-embedded Human colon cancer tissue using $\alpha\text{-Smooth}$ Muscle Actin (ACTA2) Rabbit mAb (A17910, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x. Perform high pressure antigen retrieval with 0.01 M citrate buffer (pH 6.0) prior to IF staining.