Leader in Biomolecular Solutions for Life Science



Catalog No.: A0883 10 Publications



Basic Information

Observed MW 70kDa

Calculated MW 70kDa

Category Primary antibody

Applications WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity Human, Mouse, Rat

Background

The p70/p80 autoantigen is a nuclear complex consisting of two subunits with molecular masses of approximately 70 and 80 kDa. The complex functions as a single-stranded DNA-dependent ATP-dependent helicase. The complex may be involved in the repair of nonhomologous DNA ends such as that required for double-strand break repair, transposition, and V(D)J recombination. High levels of autoantibodies to p70 and p80 have been found in some patients with systemic lupus erythematosus.

Recommended Dilutions

WB	1:3000 - 1:18000
IHC-P	1:50 - 1:200
IF/ICC	1:50 - 1:200
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID 2547

Swiss Prot P12956

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms

ML8; KU70; TLAA; CTC75; CTCBF; G22P1; Ku70

Contact


6	400-999-6126
\mathbf{X}	cn.market@abclonal.com.cn
€	www.abclonal.com.cn


Product Information

Source Rabbit **Isotype** IgG Purification Affinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

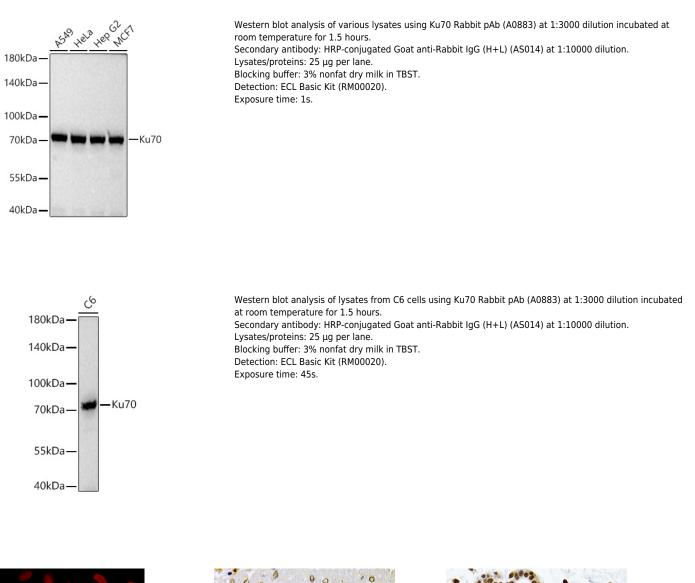


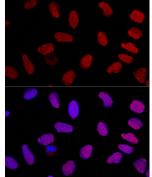


Confocal immunofluorescence analysis of U2OS cells using Ku70 Rabbit pAb (A0883) at dilution of 1:100. Blue: DAPI for nuclear staining.

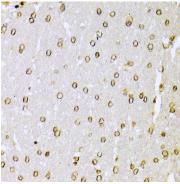

Immunohistochemistry analysis of paraffinembedded Rat brain using Ku70 Rabbit pAb (A0883) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

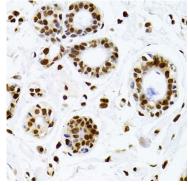

Immunohistochemistry analysis of paraffinembedded Human breast using Ku70 Rabbit pAb (A0883) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

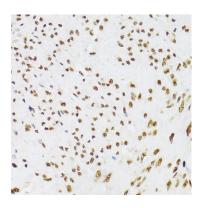

Immunohistochemistry analysis of paraffinembedded Human leiomyoma of uterus using Ku70 Rabbit pAb (A0883) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.



Immunohistochemistry analysis of paraffinembedded Mouse lung using Ku70 Rabbit pAb (A0883) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.