

Cardiac troponin T (TNNT2) Rabbit pAb

Catalog No.: A1126SP **4 Publications**

Basic Information

Observed MW

40 kDa

Calculated MW

31-36 kDa

Category

Primary antibody

Applications

WB,IF-P,IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

This gene encodes the cardiac isoform of troponin T. The encoded protein is the tropomyosin-binding subunit of the troponin complex, which is located on the thin filament of striated muscles and regulates muscle contraction in response to alterations in intracellular calcium ion concentration. Mutations in this gene have been associated with familial hypertrophic cardiomyopathy as well as with dilated cardiomyopathy.

Recommended Dilutions

WB 1:2000 - 1:6000

IF-P 1:200 - 1:2000

IHC-P 1:5000 - 1:20000

ELISA Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements. For high-ratio antibody dilutions (≥1:10000) a sequential dilution method is strongly recommended to ensure measurement accuracy.

Immunogen Information

Gene ID

7139

Swiss Prot

P45379

Immunogen

This information is considered to be commercially sensitive.

Synonyms

CMH2; RCM3; TnTC; cTnT; CMD1D; CMPD2; LVNC6; Cardiac troponin T (TNNT2)

Product Information

Source

Rabbit

Isotype

IgG

Purification

Affinity purification

Storage

Store at -20°C

Contact


☎ | 400-999-6126

✉ | cn.market@abclonal.com.cn




www.abclonal.com.cn


Validation Data

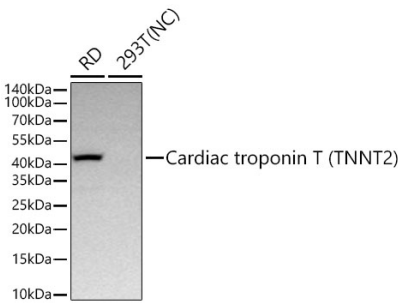


Western blot analysis of various lysates using Cardiac troponin T (TNNT2) Rabbit pAb (A1126SP) at 1:6000 dilution incubated overnight at 4°C.
Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
Lysates/proteins: 25 µg per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Negative control (NC): 293T.
Exposure time: 60 s.

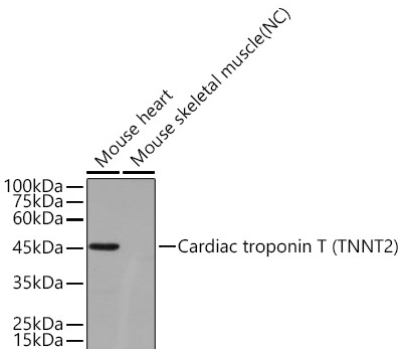


Western blot analysis of various lysates using Cardiac troponin T (TNNT2) Rabbit pAb (A1126SP) at 1:1000 dilution incubated overnight at 4°C.
Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
Lysates/proteins: 25 µg per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Negative control (NC): Mouse skeletal muscle.
Exposure time: 10 s.

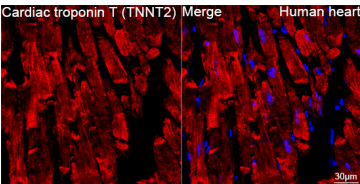

Confocal imaging of paraffin-embedded Human heart tissue using Cardiac troponin T (TNNT2) Rabbit pAb (A1126SP, dilution 1:200) followed by a further incubation with Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

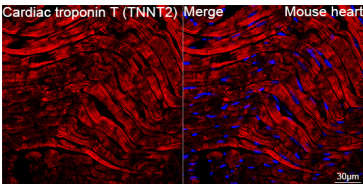


Confocal imaging of paraffin-embedded Mouse heart tissue using Cardiac troponin T (TNNT2) Rabbit pAb (A1126SP, dilution 1:200) followed by a further incubation with Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

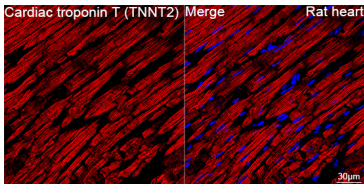


Confocal imaging of paraffin-embedded Rat heart tissue using Cardiac troponin T (TNNT2) Rabbit pAb (A1126SP, dilution 1:200) followed by a further incubation with Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

Validation Data

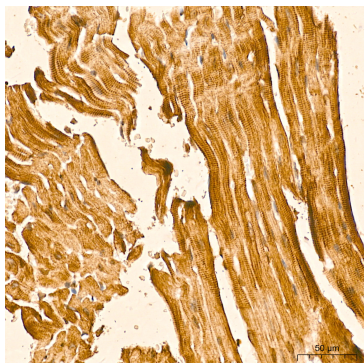

Immunohistochemistry analysis of paraffin-embedded Mouse heart tissue using Cardiac troponin T (TNNT2) Rabbit pAb (A1126SP) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

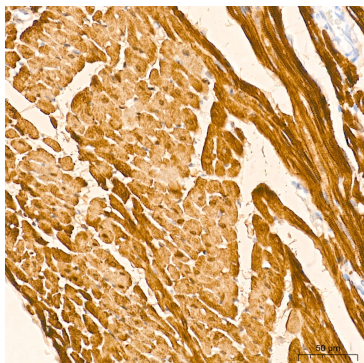

Immunohistochemistry analysis of paraffin-embedded Rat heart tissue using Cardiac troponin T (TNNT2) Rabbit pAb (A1126SP) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.