Leader in Biomolecular Solutions for Life Science

# Cytokeratin 5 Rabbit mAb

Catalog No.: A11396 Recombinant 5 Publications



## **Basic Information**

**Observed MW** 62kDa

**Calculated MW** 62kDa

Category Primary antibody

Applications WB,IHC-P,IF/ICC,ELISA

**Cross-Reactivity** Human, Mouse, Rat

**CloneNo number** ARC0585

# Background

The protein encoded by this gene is a member of the keratin gene family. The type II cytokeratins consist of basic or neutral proteins which are arranged in pairs of heterotypic keratin chains coexpressed during differentiation of simple and stratified epithelial tissues. This type II cytokeratin is specifically expressed in the basal layer of the epidermis with family member KRT14. Mutations in these genes have been associated with a complex of diseases termed epidermolysis bullosa simplex. The type II cytokeratins are clustered in a region of chromosome 12q12-q13.

## **Recommended Dilutions**

WB	1:1000 - 1:2000
IHC-P	1:400 - 1:4000
IF/ICC	1:200 - 1:400
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

## Immunogen Information

Gene ID 3852

**Swiss Prot** P13647

#### Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

#### Synonyms

K5; KRT5; DDD; DDD1; EBS1; EBS2; EBS2A; EBS2B; EBS2C; EBS2D; EBS2E; EBS2F; KRT5A; Cytokeratin 5 (CK5)

## Contact

6	400-999-6126
$\times$	cn.market@abclonal.com.cn
€	www.abclonal.com.cn

# **Product Information**


Source Rabbit


Isotype lgG


**Purification** Affinity purification


#### Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

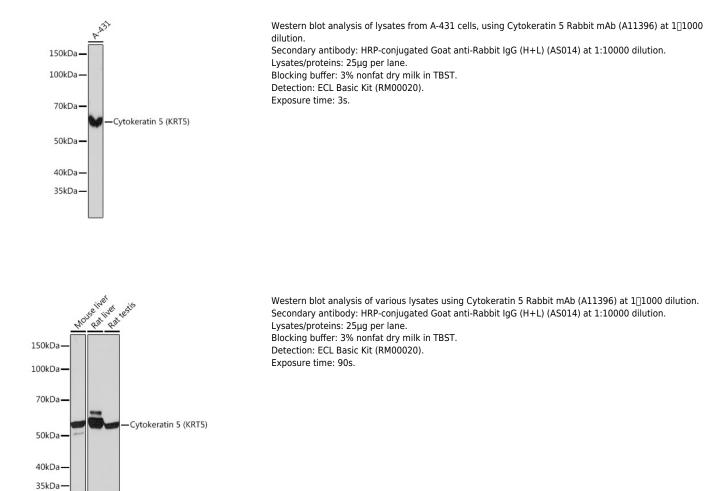

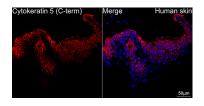

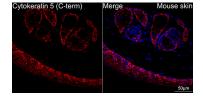


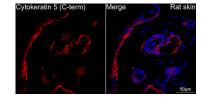


Confocal imaging of paraffin-embedded Human skin using Cytokeratin 5 Rabbit mAb (A11396, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x.Perform high pressure antigen retrieval with 0.01M citrate buffer (pH 6.0) prior to IF staining. Confocal imaging of paraffin-embedded Mouse skin using Cytokeratin 5 Rabbit mAb (A11396, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x.Perform high pressure antigen retrieval with 0.01M citrate buffer (pH 6.0) prior to IF staining. Confocal imaging of paraffin-embedded Rat skin using Cytokeratin 5 Rabbit mAb (A11396, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x.Perform high pressure antigen retrieval with 0.01M citrate buffer (pH 6.0) prior to IF staining.




Immunohistochemistry analysis of paraffinembedded Human breast tissue using Cytokeratin 5 Rabbit mAb (A11396) at a dilution of 1:2000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

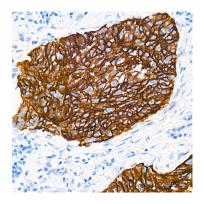


Immunohistochemistry analysis of paraffinembedded Human lung squamous carcinoma tissue using Cytokeratin 5 Rabbit mAb (A11396) at a dilution of 1:2000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

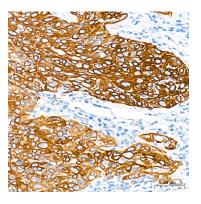

Immunohistochemistry analysis of paraffinembedded Human cervix cancer tissue using Cytokeratin 5 Rabbit mAb (A11396) at a dilution of 1:2000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

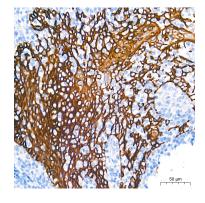

Immunohistochemistry analysis of paraffinembedded Human tonsil tissue using Cytokeratin 5 Rabbit mAb (A11396) at a dilution of 1:2000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

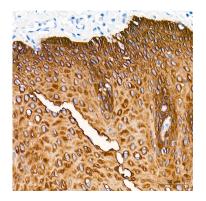

Immunohistochemistry analysis of paraffinembedded Human esophagus tissue using Cytokeratin 5 Rabbit mAb (A11396) at a dilution of 1:2000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.