NQO2 Rabbit pAb

Catalog No.: A11746



Basic Information

Observed MW

26kDa

Calculated MW

26kDa

Category

Primary antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Mouse

Background

This gene encodes a member of the thioredoxin family of enzymes. It is a cytosolic and ubiquitously expressed flavoprotein that catalyzes the two-electron reduction of quinone substrates and uses dihydronicotinamide riboside as a reducing coenzyme. Mutations in this gene have been associated with neurodegenerative diseases and several cancers. Alternative splicing results in multiple transcript variants.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:50 - 1:200

Recommended starting concentration is 1 µg/mL.

Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID Swiss Prot4835
P16083

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

QR2; DHQV; DIA6; NMOR2; NQO2

Contact

<u>a</u>		400-999-6126
\bowtie		cn.market@abclonal.com.cn
\overline{a}	Т	www.ahclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

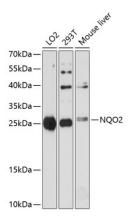


Western blot analysis of various lysates using NQO2 Rabbit pAb (A11746) at 1:3000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit lgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: $25\mu g$ per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 90s.

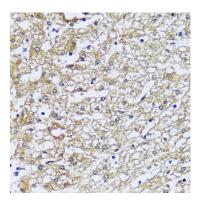

Immunohistochemistry analysis of paraffinembedded Human liver using NQO2 Rabbit pAb (A11746) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

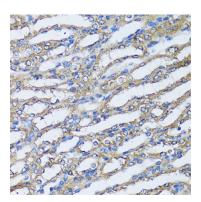

Immunohistochemistry analysis of paraffinembedded Mouse kidney using NQO2 Rabbit pAb (A11746) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.