INADL Rabbit pAb

Catalog No.: A12063 2 Publications



Basic Information

Observed MW

170kDa

Calculated MW

196kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse

Background

This gene encodes a protein with multiple PDZ domains. PDZ domains mediate protein-protein interactions, and proteins with multiple PDZ domains often organize multimeric complexes at the plasma membrane. This protein localizes to tight junctions and to the apical membrane of epithelial cells. A similar protein in Drosophila is a scaffolding protein which tethers several members of a multimeric signaling complex in photoreceptors.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:50 - 1:200

Immunogen Information

Gene IDSwiss Prot
10207
Q8NI35

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 1-100 of human INADL ($NP_795352.2$).

Synonyms

Cipp; INADL; hINADL; InaD-like

Contact

<u>a</u>		400-999-6126
\bowtie		cn.market@abclonal.com.cn
$\overline{\triangle}$	ī	www.ahclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20 $^{\circ}\text{C}.$ Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

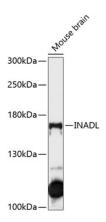


Western blot analysis of lysates from mouse brain, using INADL Rabbit pAb (A12063) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 90s.

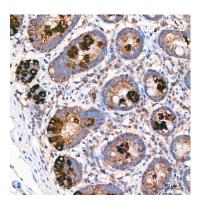

Immunohistochemistry analysis of INADL in paraffin-embedded Human colon tissue using INADL Rabbit pAb (A12063) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer(pH 6.0) prior to IHC staining.

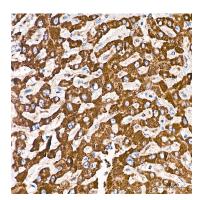

Immunohistochemistry analysis of INADL in paraffin-embedded Human liver tissue using INADL Rabbit pAb (A12063) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer(pH 6.0) prior to IHC staining.