# VCAM1 Rabbit mAb

Catalog No.: A19131 Recombinant 7 Publications



## **Basic Information**

### **Observed MW**

110kDa

### **Calculated MW**

81kDa

### Category

Primary antibody

### **Applications**

WB,IF/ICC,ELISA

### **Cross-Reactivity**

Human, Mouse, Rat

#### CloneNo number

ARC0312

## **Background**

This gene is a member of the Ig superfamily and encodes a cell surface sialoglycoprotein expressed by cytokine-activated endothelium. This type I membrane protein mediates leukocyte-endothelial cell adhesion and signal transduction, and may play a role in the development of artherosclerosis and rheumatoid arthritis. Three alternatively spliced transcripts encoding different isoforms have been described for this gene.

## **Recommended Dilutions**

WB 1:1000 - 1:2000

**IF/ICC** 1:100 - 1:800

**ELISA** Recommended starting concentration is 1 μg/mL.

Please optimize the concentration based on your specific assay requirements.

## **Immunogen Information**

**Gene ID**7412

Swiss Prot
P19320

### **Immunogen**

Synthetic peptide. This information is considered to be commercially sensitive.

### **Synonyms**

CD106; INCAM-100; VCAM1

## **Contact**

<b>a</b>	400-999-6126
<b>×</b>	cn.market@abclonal.com.cn
$\overline{\Box}$	www.ahclonal.com.cn


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

#### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

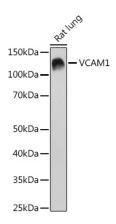


Western blot analysis of lysates from Rat lung, using VCAM1 Rabbit mAb (A19131) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit lgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins:  $25\mu g$  per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 3min.

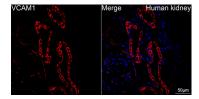

Confocal imaging of human kidney using VCAM1 Rabbit mAb (A19131,at dilution of 1:100) (Red). DAPI was used for nuclear staining (blue). Objective: 40x. Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IF staining protocol.