eNOS Rabbit mAb

Catalog No.: A20985 Recombinant 10 Publications



Basic Information

Observed MW

140kDa

Calculated MW

133kDa

Category

Primary antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC51618

Background

Nitric oxide is a reactive free radical which acts as a biologic mediator in several processes, including neurotransmission and antimicrobial and antitumoral activities. Nitric oxide is synthesized from L-arginine by nitric oxide synthases. Variations in this gene are associated with susceptibility to coronary spasm. Alternative splicing and the use of alternative promoters results in multiple transcript variants.

Recommended Dilutions

WB 1:1000 - 1:2000

IHC-P 1:200 - 1:2000

ELISA Recommended starting concentration is 1 μg/mL.

Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID4846

Swiss Prot
P29474

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

eNOS; ECNOS

Contact

a	400-999-6126
×	cn.market@abclonal.com.cn
	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

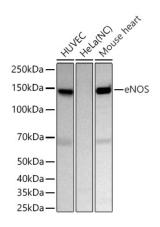

Western blot analysis of various lysates using eNOS Rabbit mAb (A20985) at 1:1000 dilution incubated at room temperature for 1.5 hours.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25 µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Negative control (NC): HeLa. Exposure time: 45 s.

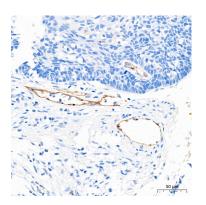


Immunohistochemistry analysis of paraffinembedded Human cervix cancer tissue using eNOS Rabbit mAb (A20985) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

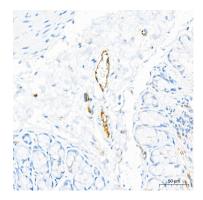


Immunohistochemistry analysis of paraffinembedded Rat colon tissue using eNOS Rabbit mAb (A20985) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

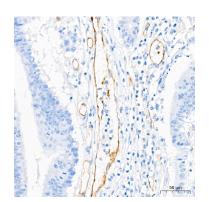

Immunohistochemistry analysis of paraffinembedded Human colon carcinoma tissue using eNOS Rabbit mAb (A20985) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

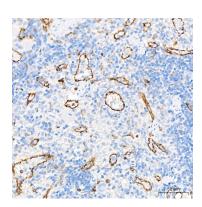

Immunohistochemistry analysis of paraffinembedded Rat spleen tissue using eNOS Rabbit mAb (A20985) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

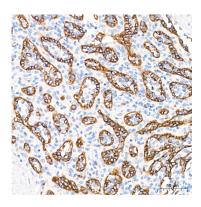

Immunohistochemistry analysis of paraffinembedded Human spleen tissue using eNOS Rabbit mAb (A20985) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.