# **GLUD1/2 Rabbit mAb**

Catalog No.: A21792 Recombinant



# **Basic Information**

### **Observed MW**

52kDa

### **Calculated MW**

61kDa

### Category

Primary antibody

### **Applications**

WB,IHC-P,IF/ICC,ELISA

#### **Cross-Reactivity**

Human, Mouse, Rat

#### CloneNo number

ARC53980

# **Background**

The protein encoded by this gene is localized to the mitochondrion and acts as a homohexamer to recycle glutamate during neurotransmission. The encoded enzyme catalyzes the reversible oxidative deamination of glutamate to alpha-ketoglutarate. This gene is intronless.

# **Recommended Dilutions**

WB 1:10000-1:130000

**IHC-P** 1:50 - 1:200

**IF/ICC** 1:50 - 1:200

**ELISA** Recommended starting

concentration is 1 µg/mL.

Please optimize the
concentration based on
your specific assay
requirements.

# **Immunogen Information**

 Gene ID
 Swiss Prot

 2746/2747
 P00367/P49448

#### **Immunogen**

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

### **Synonyms**

GLUD2; GDH2; GLUDP1; GLUD1/2

## **Contact**

	400-999-6126
<b>×</b>	cn.market@abclonal.com.cn
	www.abclonal.com.cn


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

### **Storage**

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

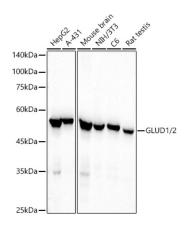


Western blot analysis of various lysates, using GLUD1/2 Rabbit mAb (A21792) at1:110000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins:  $25\mu g$  per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.

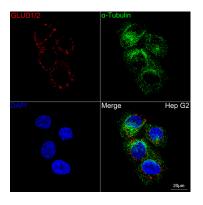


GLUD1/2 Merge Human kidney




Confocal imaging of Hep G2 cells using GLUD1/2 Rabbit mAb (A21792, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with  $\alpha$ -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

Confocal imaging of paraffin-embedded Human kidney tissue using GLUD1/2 Rabbit mAb (A21792, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x. Perform high pressure antigen retrieval with 0.01 M citrate buffer (pH 6.0) prior to IF staining.


Immunohistochemistry analysis of paraffinembedded Human colon carcinoma tissue using GLUD1/2 Rabbit mAb (A21792) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

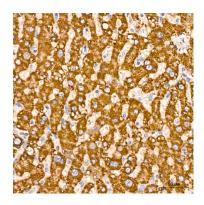


Immunohistochemistry analysis of paraffinembedded Human liver tissue using GLUD1/2 Rabbit mAb (A21792) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

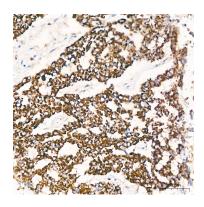


Immunohistochemistry analysis of paraffinembedded Human lung cancer tissue using GLUD1/2 Rabbit mAb (A21792) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

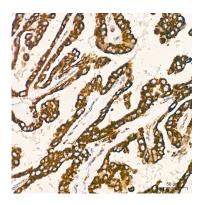

Immunohistochemistry analysis of paraffinembedded Human thyroid cancer tissue using GLUD1/2 Rabbit mAb (A21792) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

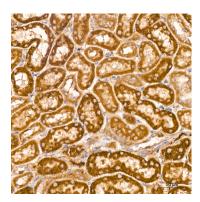

Immunohistochemistry analysis of paraffinembedded Mouse kidney tissue using GLUD1/2 Rabbit mAb (A21792) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

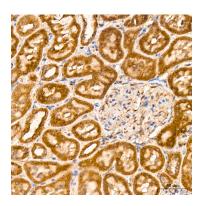

Immunohistochemistry analysis of paraffinembedded Rat kidney tissue using GLUD1/2 Rabbit mAb (A21792) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.