# TdT Rabbit mAb

Catalog No.: A22618 Recombinant



### **Basic Information**

### **Observed MW**

59kDa

### **Calculated MW**

59kDa

## Category

Primary antibody

### **Applications**

ELISA,WB,IHC-P

### **Cross-Reactivity**

Human, Mouse, Rat

#### CloneNo number

ARC55687

# **Background**

This gene is a member of the DNA polymerase type-X family and encodes a template-independent DNA polymerase that catalyzes the addition of deoxynucleotides to the 3'-hydroxyl terminus of oligonucleotide primers. In vivo, the encoded protein is expressed in a restricted population of normal and malignant pre-B and pre-T lymphocytes during early differentiation, where it generates antigen receptor diversity by synthesizing non-germ line elements (N-regions) at the junctions of rearranged Ig heavy chain and T cell receptor gene segments. Alternatively spliced transcript variants encoding different isoforms of this gene have been described.

# **Recommended Dilutions**

WB 1:10000 - 1:70000

**IHC-P** 1:50 - 1:200

# **Immunogen Information**

**Gene ID**1791

Swiss Prot
P04053

### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 10-124 of human TdT. (NP\_004079.3).

# **Synonyms**

TDT; TdT

## **Contact**

<b>a</b>	400-999-6126
<b>×</b>	cn.market@abclonal.com.cn
$\overline{\Box}$	www abclonal com cn


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

### **Storage**

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,0.05% BSA,50% glycerol,pH7.3.

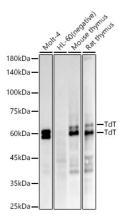


Western blot analysis of various lysates, using TdT Rabbit mAb (A22618) at 1:60000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit (RM00021).

Exposure time: 180s.

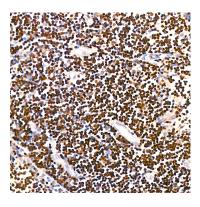

Immunohistochemistry analysis of paraffinembedded Human thymoma tissue using TdT Rabbit mAb (A22618) at a dilution of 1:200 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.