Leader in Biomolecular Solutions for Life Science

# TSG101/VPS23 Rabbit mAb

Catalog No.: A22762 Recombinant 1 Publications



### **Basic Information**

**Observed MW** Refer to figures

**Calculated MW** 44kDa

Category Primary antibody

Applications ELISA, IHC-P

**Cross-Reactivity** Mouse

**CloneNo number** ARC57589

**Recommended Dilutions** 

1:50 - 1:200

### Background

The protein encoded by this gene belongs to a group of apparently inactive homologs of ubiquitin-conjugating enzymes. The gene product contains a coiled-coil domain that interacts with stathmin, a cytosolic phosphoprotein implicated in tumorigenesis. The protein may play a role in cell growth and differentiation and act as a negative growth regulator. In vitro steady-state expression of this tumor susceptibility gene appears to be important for maintenance of genomic stability and cell cycle regulation. Mutations and alternative splicing in this gene occur in high frequency in breast cancer and suggest that defects occur during breast cancer tumorigenesis and/or progression.

## **Immunogen Information**

IHC-P

#### Gene ID 7251

**Swiss Prot** Q99816

#### Immunogen

A synthesized peptide derived from human TSG101/VPS23.

#### **Synonyms**

TSG10; VPS23; TSG101/VPS23

#### 400-999-6126 8 cn.market@abclonal.com.cn $\sim$ Ð www.abclonal.com.cn

### **Product Information**


Source Rabbit

Isotype lgG


Purification Affinity purification

#### Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.05% proclin300,0.05% BSA,50% glycerol,pH7.3.

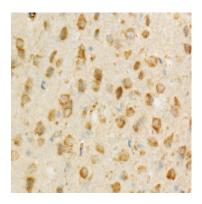

Immunohistochemistry analysis of paraffinembedded mouse brain using TSG101/VPS23 Rabbit mAb (A22762) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

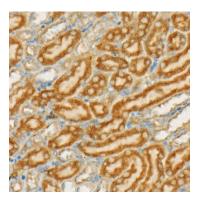

Immunohistochemistry analysis of paraffinembedded mouse kidney using TSG101/VPS23 Rabbit mAb (A22762) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.