CPS1 Rabbit mAb

Catalog No.: A24810 Recombinant



Basic Information

Observed MW

165kDa

Calculated MW

165kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC63814

Background

The mitochondrial enzyme encoded by this gene catalyzes synthesis of carbamoyl phosphate from ammonia and bicarbonate. This reaction is the first committed step of the urea cycle, which is important in the removal of excess urea from cells. The encoded protein may also represent a core mitochondrial nucleoid protein. Three transcript variants encoding different isoforms have been found for this gene. The shortest isoform may not be localized to the mitochondrion. Mutations in this gene have been associated with carbamoyl phosphate synthetase deficiency, susceptibility to persistent pulmonary hypertension, and susceptibility to venoocclusive disease after bone marrow transplantation.

Recommended Dilutions

WB 1:2000 - 1:5000

IHC-P 1:200 - 1:500

Immunogen Information

Gene IDSwiss Prot
1373
P31327

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1356-1482 of human CPS1 (NP_001116105.1).

Synonyms

PHN; GATD6; CPSASE1; CPS1

Contact

6		400-999-6126
\bowtie		cn.market@abclonal.com.cn
•	T	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,0.05% BSA,50% glycerol,pH7.3.

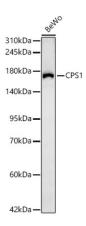


Western blot analysis of lysates from BeWo cells, using CPS1 Rabbit mAb (A24810) at 1:4000 dilution. Secondary antibody: HRP Goat Anti-Rabbit lgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25ug per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 60s.

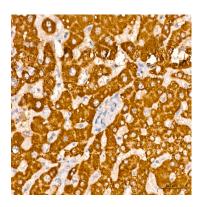

Immunohistochemistry analysis of CPS1 in paraffin-embedded human liver tissue using CPS1 Rabbit mAb (A24810) at dilution of 1:300 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.