[KO Validated] MLH1 Rabbit mAb

ABclonal www.abclonal.com

Catalog No.: A25641 KO Validated Recombinant

1 Publications

Basic Information

Observed MW

85kDa/85 kDa/85 kDa

Calculated MW

85kDa

Category

Primary antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Mouse

CloneNo number

ARC53543

Background

The protein encoded by this gene can heterodimerize with mismatch repair endonuclease PMS2 to form MutL alpha, part of the DNA mismatch repair system. When MutL alpha is bound by MutS beta and some accessory proteins, the PMS2 subunit of MutL alpha introduces a single-strand break near DNA mismatches, providing an entry point for exonuclease degradation. The encoded protein is also involved in DNA damage signaling and can heterodimerize with DNA mismatch repair protein MLH3 to form MutL gamma, which is involved in meiosis. This gene was identified as a locus frequently mutated in hereditary nonpolyposis colon cancer (HNPCC).

Recommended Dilutions

WB 1:8000 - 1:16000

1:100 - 1:500 **IHC-P**

ELISA Recommended starting concentration is 1 µg/mL.

Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID Swiss Prot 4292 P40692

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms

FCC2; COCA2; HNPCC; MLH-1; hMLH1; HNPCC2; LYNCH2; MMRCS1

Contact

a	400-999-6126
×	cn.market@abclonal.com.cn
\overline{a}	www.abclonal.com.cn


Product Information

Source Isotype **Purification** Rabbit IgG Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

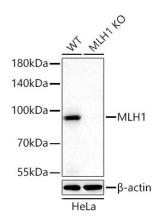


Western blot analysis of lysates from wild type (WT) and MLH1 knockout (KO) HeLa cells using [KO Validated] MLH1 Rabbit mAb (A25641) at 1:8000 dilution incubated at room temperature for 1.5 hours. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25 µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 60 s.

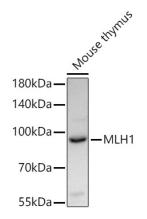

Western blot analysis of lysates from Mouse thymus using [KO Validated] MLH1 Rabbit mAb (A25641) at 1:8000 dilution incubated at room temperature for 1.5 hours.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25 µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 20 s.

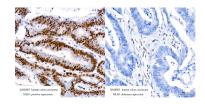

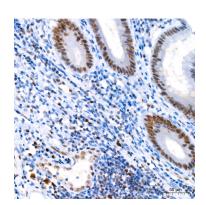


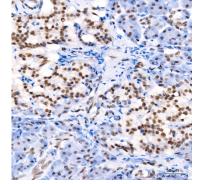

Immunohistochemistry analysis of paraffinembedded Human colon carcinoma tissue using [KO Validated] MLH1 Rabbit mAb (A25641) at a dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

Immunohistochemistry analysis of paraffinembedded Human appendix tissue using [KO Validated] MLH1 Rabbit mAb (A25641) at a dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

Immunohistochemistry analysis of paraffinembedded Human pancreas tissue using [KO Validated] MLH1 Rabbit mAb (A25641) at a dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

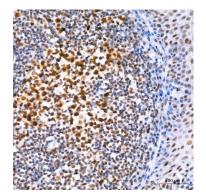







Immunohistochemistry analysis of paraffinembedded Human tonsil tissue using [KO Validated] MLH1 Rabbit mAb (A25641) at a dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

Immunohistochemistry analysis of paraffinembedded Human HeLa and HeLa-MLH1-KO cell line using [KO Validated] MLH1 Rabbit mAb (A25641) at a dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

Immunohistochemistry analysis of paraffinembedded Human colon carcinoma using [KO Validated] MLH1 Rabbit mAb (A4858) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.