GLDC Rabbit PolymAb®

Catalog No.: A26427PM



Basic Information

Observed MW

113kDa

Calculated MW

113kDa

Category

Primary antibody

Applications

WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

Degradation of glycine is brought about by the glycine cleavage system, which is composed of four mitochondrial protein components: P protein (a pyridoxal phosphate-dependent glycine decarboxylase), H protein (a lipoic acid-containing protein), T protein (a tetrahydrofolate-requiring enzyme), and L protein (a lipoamide dehydrogenase). The protein encoded by this gene is the P protein, which binds to glycine and enables the methylamine group from glycine to be transferred to the T protein. Defects in this gene are a cause of nonketotic hyperglycinemia (NKH).

Recommended Dilutions

WB 1:25000 - 1:150000

IHC-P 1:4000 - 1:16000

IF/ICC 1:2000 - 1:8000

ELISA Recommended starting

concentration is 1 µg/mL.

Please optimize the
concentration based on
your specific assay
requirements.

Immunogen Information

Gene ID2731

Swiss Prot
P23378

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

GCE; GCSP; HYGN1

Contact

a		400-999-6126
\bowtie		cn.market@abclonal.com.cn
\odot	T	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

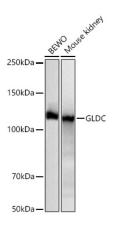

Western blot analysis of various lysates using GLDC Rabbit PolymAb® (A26427PM) at 1:25000 dilution incubated overnight at 4° C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25 µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 1s.

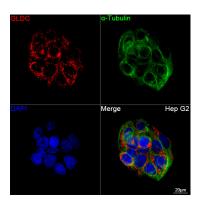


Confocal imaging of Hep G2 cells using GLDC Rabbit PolymAb® (A26427PM, dilution 1:2000) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

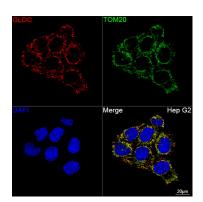


Confocal imaging of Hep G2 cells using GLDC Rabbit PolymAb® (A26427PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with TOMM20 Mouse mAb (A27799, dilution 1:200) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution

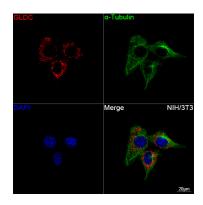


Confocal imaging of NIH/3T3 cells using GLDC Rabbit PolymAb® (A26427PM, dilution 1:2000) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

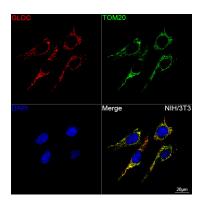

Confocal imaging of NIH/3T3 cells using GLDC Rabbit PolymAb® (A26427PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with TOMM20 Mouse mAb (A27799, dilution 1:200) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution

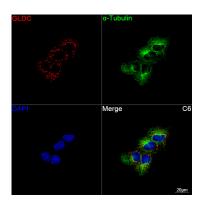


Confocal imaging of C6 cells using GLDC Rabbit PolymAb® (A26427PM, dilution 1:2000) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

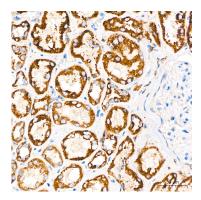

Immunohistochemistry analysis of paraffinembedded Human kidney tissue using GLDC Rabbit PolymAb® (A26427PM) at a dilution of 1:4000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

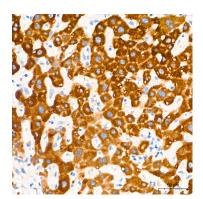


Immunohistochemistry analysis of paraffinembedded Human liver tissue using GLDC Rabbit PolymAb® (A26427PM) at a dilution of 1:4000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.



Immunohistochemistry analysis of paraffinembedded Mouse liver tissue using GLDC Rabbit PolymAb® (A26427PM) at a dilution of 1:4000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

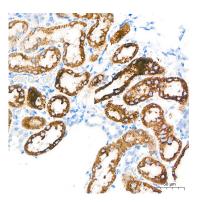

Immunohistochemistry analysis of paraffinembedded Rat kidney tissue using GLDC Rabbit PolymAb® (A26427PM) at a dilution of 1:4000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.