Fibrillarin Rabbit mAb

Catalog No.: A27818 Recombinant



Basic Information

Observed MW

37KD

Calculated MW

34kDa

Category

Primary antibody

Applications

WB,IP,IF/ICC,IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC71609

Background

This gene product is a component of a nucleolar small nuclear ribonucleoprotein (snRNP) particle thought to participate in the first step in processing preribosomal RNA. It is associated with the U3, U8, and U13 small nuclear RNAs and is located in the dense fibrillar component (DFC) of the nucleolus. The encoded protein contains an N-terminal repetitive domain that is rich in glycine and arginine residues, like fibrillarins in other species. Its central region resembles an RNA-binding domain and contains an RNP consensus sequence. Antisera from approximately 8% of humans with the autoimmune disease scleroderma recognize fibrillarin.

Recommended Dilutions

WB 1:25000 - 1:100000

IP 0.5μg-4μg antibody for 200μg-400μg extracts of

whole cells

IF/ICC 1:200 - 1:400

IHC-P 1:1000 - 1:4000

ELISA Recommended starting

concentration is 1 µg/mL.

Please optimize the
concentration based on
your specific assay
requirements.

Contact

<u>a</u>	400-999-6126
\bowtie	cn.market@abclonal.com.cn
•	www.abclonal.com.cn

Immunogen Information

Gene ID Swiss Prot

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms

FIB; FLRN; Nop1; RNU3IP1


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

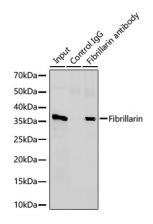

Immunoprecipitation of Fibrillarin from 300 μg extracts of 293T cells was performed using 0.5 μg of Fibrillarin Rabbit mAb (A27818). Rabbit IgG isotype control (AC005) was used to precipitate the Control IgG sample. IP samples were eluted with 1X reducing Laemmli Buffer. The Input lane represents 10% of the total input. Western blot analysis of immunoprecipitates was conducted using Fibrillarin Rabbit mAb (A27818) at a dilution of 1:5000.

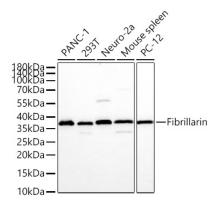


Western blot analysis of various lysates using Fibrillarin Rabbit mAb (A27818) at 1:25000 dilution incubated overnight at 4° C.

Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: 25 μ g per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 20s.

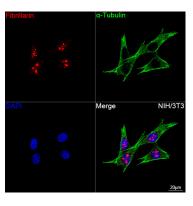

Confocal imaging of NIH/3T3 cells using Fibrillarin Rabbit mAb (A27818, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

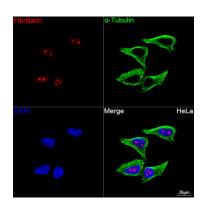

Confocal imaging of HeLa cells using Fibrillarin Rabbit mAb (A27818, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

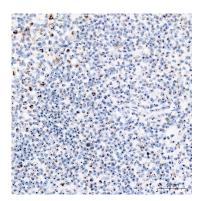


Immunohistochemistry analysis of paraffinembedded Human tonsil tissue using Fibrillarin Rabbit mAb (A27818) at a dilution of 1:3000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

Validation Data



Immunohistochemistry analysis of paraffinembedded Mouse liver tissue using Fibrillarin Rabbit mAb (A27818) at a dilution of 1:3000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

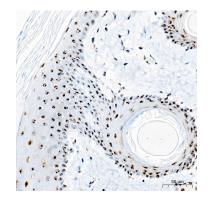

Immunohistochemistry analysis of paraffinembedded Rat skin tissue using Fibrillarin Rabbit mAb (A27818) at a dilution of 1:3000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.