


Validation Data

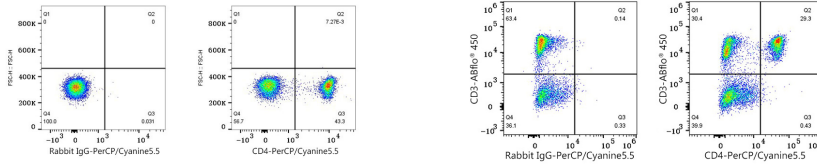

Flow cytometry: 1×10^6 Human PBMC were surface-stained with PerCP/Cyanine5.5 Rabbit IgG isotype control (A25620, 5 μ l/Test, left) or PerCP/Cyanine5.5 Rabbit anti-Human/Monkey CD4 mAb (A27881, 5 μ l/Test, right). Cells in the lymphocyte gate were used for analysis.

Flow cytometry: 1×10^6 Cynomolgus peripheral blood mononuclear cells were surface-stained with ABflo® 450 Rabbit anti-Human/Monkey CD3 mAb (A27177, 5 μ l/Test) and PerCP/Cyanine5.5 Rabbit IgG isotype control (A25620, 5 μ l/Test, left) or PerCP/Cyanine5.5 Rabbit anti-Human/Monkey CD4 mAb (A27881, 5 μ l/Test, right). Cells in the lymphocyte gate were used for analysis.