



Validation Data

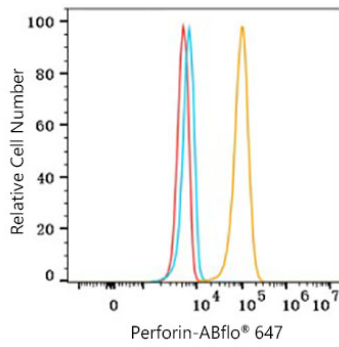

Flow cytometry: 1×10^6 KARPAS-299 cells were intracellularly-stained with ABflo® 647 Rabbit anti-Human Perforin mAb (A28172, 5 μ l/Test, orange line) or ABflo® 647 Rabbit IgG isotype control (A22070, 5 μ l/Test, blue line). Non-fluorescently stained cells were used as blank control (red line).

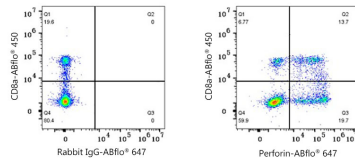

Flow cytometry: 1×10^6 Human PBMC were surface-stained with ABflo® 450 Rabbit anti-Human/Monkey CD8a mAb (A27220, 5 μ l/Test) and then intracellularly-stained with ABflo® 647 Rabbit IgG isotype control (A22070, 5 μ l/Test, left) or ABflo® 647 Rabbit anti-Human Perforin mAb (A28172, 5 μ l/Test, Right). Cells in the lymphocyte gate were used for analysis.