APC Rabbit pAb

Catalog No.: A2818 2 Publications



Basic Information

Observed MW

150kDa

Calculated MW

312kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P,IF/ICC

Cross-Reactivity

Human, Mouse, Rat

Background

This gene encodes a tumor suppressor protein that acts as an antagonist of the Wnt signaling pathway. It is also involved in other processes including cell migration and adhesion, transcriptional activation, and apoptosis. Defects in this gene cause familial adenomatous polyposis (FAP), an autosomal dominant pre-malignant disease that usually progresses to malignancy. Mutations in the APC gene have been found to occur in most colorectal cancers, where disease-associated mutations tend to be clustered in a small region designated the mutation cluster region (MCR) and result in a truncated protein product.

Recommended Dilutions

WB	1:500 - 1:1000
IHC-P	1:50 - 1:200
IF/ICC	1:50 - 1:200

Immunogen Information

Gene ID	Swiss Prot	
324	P25054	

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 2600-2760 of human APC (NP_000029.2).

Synonyms

GS; DP2; DP3; BTPS2; DESMD; DP2.5; PPP1R46; APC

Contact

a		400-999-6126
\bowtie		cn.market@abclonal.com.cn
•	T	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

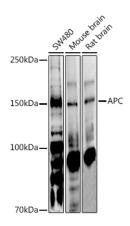


Western blot analysis of extracts of various cell lines, using APC antibody (A2818) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 180s.

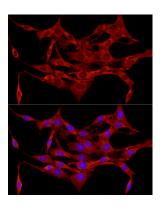

Immunofluorescence analysis of PC-12 cells using APC antibody (A2818) at dilution of 1:100. Blue: DAPI for nuclear staining.

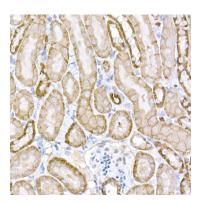

Immunohistochemistry analysis of paraffinembedded mouse kidney using APC Rabbit pAb (A2818) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.