SNAP25 Rabbit mAb

Catalog No.: A28343 Recombinant



Basic Information

Observed MW

25 kDa

Calculated MW

23 kDa

Category

Primary antibody

Applications

WB,IF-P,IHC-P,ELISA

Cross-Reactivity

Mouse, Rat

CloneNo number

ARC75855

Background

Synaptic vesicle membrane docking and fusion is mediated by SNAREs (soluble Nethylmaleimide-sensitive factor attachment protein receptors) located on the vesicle membrane (v-SNAREs) and the target membrane (t-SNAREs). The assembled v-SNARE/t-SNARE complex consists of a bundle of four helices, one of which is supplied by v-SNARE and the other three by t-SNARE. For t-SNAREs on the plasma membrane, the protein syntaxin supplies one helix and the protein encoded by this gene contributes the other two. Therefore, this gene product is a presynaptic plasma membrane protein involved in the regulation of neurotransmitter release. Two alternative transcript variants encoding different protein isoforms have been described for this gene.

Recommended Dilutions

WB 1:2000 - 1:10000

IF-P 1:200 - 1:1000

IHC-P 1:2000 - 1:8000

ELISA Recommended starting

concentration is 1 µg/mL.

Please optimize the concentration based on your specific assay requirements. For highratio antibody dilutions (≥1:10000)□a sequential dilution method is strongly recommended to ensure measurement

accuracy.

Immunogen Information

Gene IDSwiss Prot
6616
P60880

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

SUP; RIC4; SEC9; SNAP; CMS18; RIC-4; SNAP-25; bA416N4.2; dJ1068F16.2

Product Information

SourceIsotypePurificationRabbitIgGAffinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

Contact

2	400-999-6126
\bowtie	cn.market@abclonal.com.cr
•	www.abclonal.com.cr

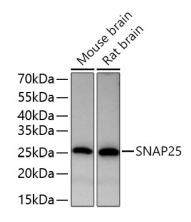

Western blot analysis of various lysates using SNAP25 Rabbit mAb (A28343) at 1:5000 dilution incubated overnight at 4° C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25 µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 45 s.

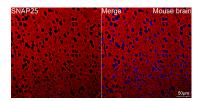

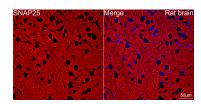


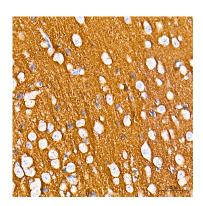


Confocal imaging of paraffin-embedded Mouse brain tissue using SNAP25 Rabbit mAb (A28343, dilution 1:200) followed by a further incubation with Cy3-conjugated Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

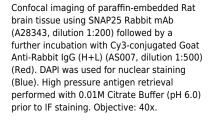


Immunohistochemistry analysis of paraffinembedded Mouse brain tissue using SNAP25 Rabbit mAb (A28343) at a dilution of 1:5000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

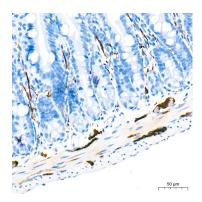

Immunohistochemistry analysis of paraffinembedded Mouse intestine tissue using SNAP25 Rabbit mAb (A28343) at a dilution of 1:5000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

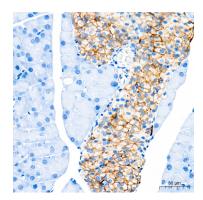


Immunohistochemistry analysis of paraffinembedded Mouse pancreas tissue using SNAP25 Rabbit mAb (A28343) at a dilution of 1:5000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

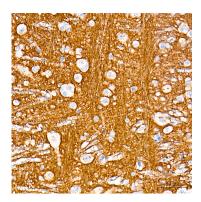


Immunohistochemistry analysis of paraffinembedded Rat brain tissue using SNAP25 Rabbit mAb (A28343) at a dilution of 1:5000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

Validation Data

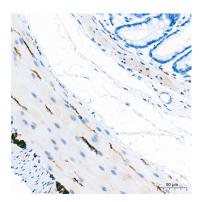

Immunohistochemistry analysis of paraffinembedded Rat colon tissue using SNAP25 Rabbit mAb (A28343) at a dilution of 1:5000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

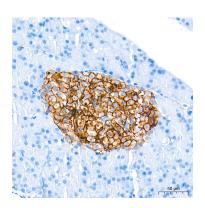

Immunohistochemistry analysis of paraffinembedded Rat pancreas tissue using SNAP25 Rabbit mAb (A28343) at a dilution of 1:5000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.