

EHMT2/G9A Rabbit mAb

Catalog No.: A29495 **Recombinant**

Basic Information

Observed MW

160 kDa/170 kDa

Calculated MW

19 kDa/129 kDa/132 kDa

Category

Primary antibody

Applications

WB,IP,IF/ICC,IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC3375

Background

This gene encodes a methyltransferase that methylates lysine residues of histone H3. Methylation of H3 at lysine 9 by this protein results in recruitment of additional epigenetic regulators and repression of transcription. This gene was initially thought to be two different genes, NG36 and G9a, adjacent to each other in the HLA locus. Alternative splicing results in multiple transcript variants.

Recommended Dilutions

WB 1:2000 - 1:20000

IP 0.5 µg - 4 µg antibody for
400 µg - 600 µg extracts of
whole cells

IF/ICC 1:1000 - 1:4000

IHC-P 1:200 - 1:1000

ELISA Recommended starting
concentration is 1 µg/mL.
Please optimize the
concentration based on your
specific assay requirements.
For high-ratio antibody
dilutions ($\geq 1:10000$) a
sequential dilution method
is strongly recommended to
ensure measurement
accuracy.

Immunogen Information

Gene ID

10919

Swiss Prot

Q96KQ7

Immunogen

This information is considered to be commercially sensitive.

Synonyms

G9A; BAT8; GAT8; NG36; KMT1C; C6orf30

Product Information

Source

Rabbit

Isotype

IgG

Purification

Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide, 0.05% BSA, 50% glycerol, pH7.3.


Contact

☎ | 400-999-6126


✉ | cn.market@abclonal.com.cn

🌐 | www.abclonal.com.cn


Validation Data

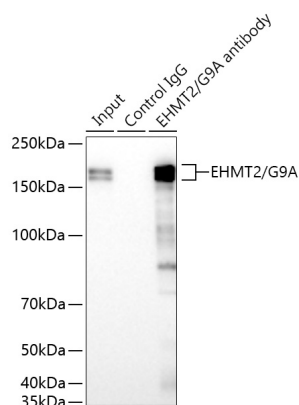


Immunoprecipitation of EHMT2/G9A from 600 µg extracts of Jurkat cells was performed using 1 µg of EHMT2/G9A Rabbit mAb (A29495). Rabbit Control IgG (AC005) was used to precipitate the Control IgG sample. IP samples were eluted with 1x Laemmli Buffer. The Input lane represents 10% of the total input. Western blot analysis of immunoprecipitates was conducted using EHMT2/G9A Rabbit mAb (A29495) at a dilution of 1:5000.

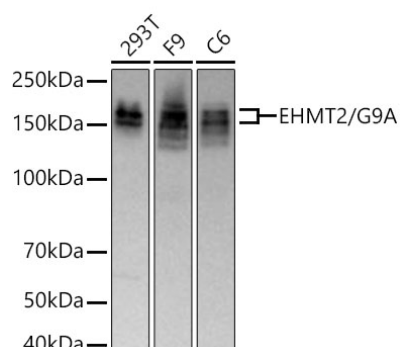


Western blot analysis of various lysates using EHMT2/G9A Rabbit mAb (A29495) at 1:20000 dilution incubated overnight at 4°C. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: 25 µg per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020). Exposure time: 90 s.

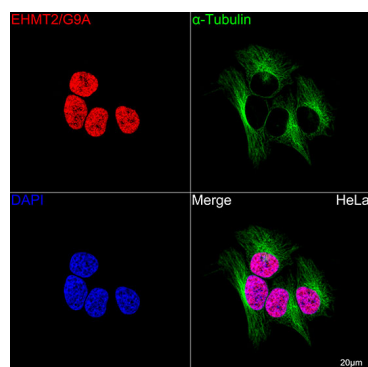

Confocal imaging of HeLa cells using EHMT2/G9A Rabbit mAb (A29495, dilution 1:2000) followed by a further incubation with Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α-Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) (AS037, dilution 1:800) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

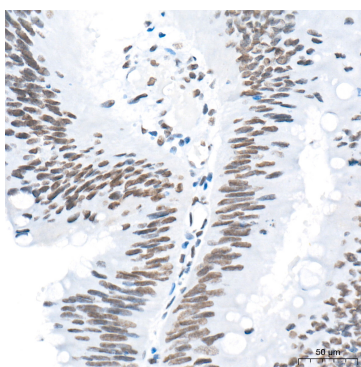


Immunohistochemistry analysis of paraffin-embedded Human colon carcinoma tissue using EHMT2/G9A Rabbit mAb (A29495) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

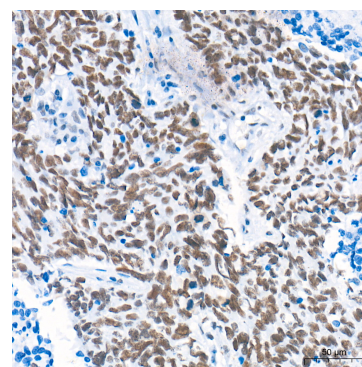


Immunohistochemistry analysis of paraffin-embedded Human lung cancer tissue using EHMT2/G9A Rabbit mAb (A29495) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.


Validation Data

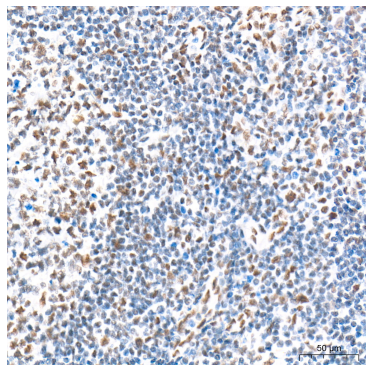

Immunohistochemistry analysis of paraffin-embedded Human tonsil tissue using EHMT2/G9A Rabbit mAb (A29495) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

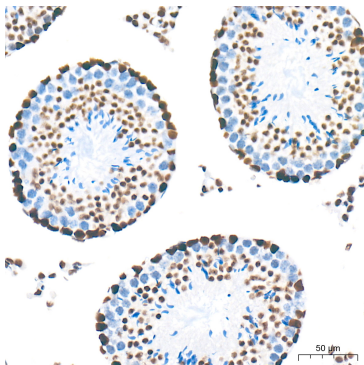

Immunohistochemistry analysis of paraffin-embedded Mouse testis tissue using EHMT2/G9A Rabbit mAb (A29495) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

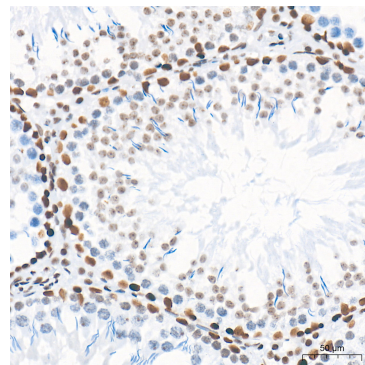

Immunohistochemistry analysis of paraffin-embedded Rat testis tissue using EHMT2/G9A Rabbit mAb (A29495) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.