FUT6 Rabbit pAb

Catalog No.: A3868



Basic Information

Observed MW

42kDa

Calculated MW

42kDa

Category

Primary antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human

Background

The protein encoded by this gene is a Golgi stack membrane protein that is involved in the creation of sialyl-Lewis X, an E-selectin ligand. Mutations in this gene are a cause of fucosyltransferase-6 deficiency. Two transcript variants encoding the same protein have been found for this gene.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:100 - 1:200

ELISA Recommended starting concentration is 1 μg/mL.

Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID2528

Swiss Prot
P51993

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

FT1A; FCT3A; Fuc-TVI; FucT-VI; FUT6

Contact

| a | | 400-999-6126 |
|-----------|---|---------------------------|
| \bowtie | | cn.market@abclonal.com.cn |
| • | ī | www.abclonal.com.cn |


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

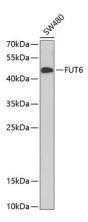


Western blot analysis of lysates from SW480 cells, using FUT6 Rabbit pAb (A3868) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 90s.

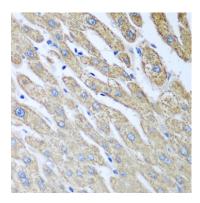

Immunohistochemistry analysis of paraffinembedded Human liver using FUT6 Rabbit pAb (A3868) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with $0.01 M\ PBS\ Buffer\ (pH\ 7.2)\ prior\ to\ IHC$ staining.

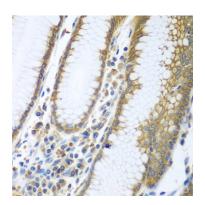

Immunohistochemistry analysis of paraffinembedded Human stomach using FUT6 Rabbit pAb (A3868) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.