# **ACTC1 Rabbit pAb**

Catalog No.: A6507



## **Basic Information**

#### **Observed MW**

42kDa

#### **Calculated MW**

42kDa

### Category

Primary antibody

### **Applications**

ELISA,WB,IHC-P

### **Cross-Reactivity**

Human, Mouse, Rat

# **Background**

Actins are highly conserved proteins that are involved in various types of cell motility. Polymerization of globular actin (G-actin) leads to a structural filament (F-actin) in the form of a two-stranded helix. Each actin can bind to four others. The protein encoded by this gene belongs to the actin family which is comprised of three main groups of actin isoforms, alpha, beta, and gamma. The alpha actins are found in muscle tissues and are a major constituent of the contractile apparatus. Defects in this gene have been associated with idiopathic dilated cardiomyopathy (IDC) and familial hypertrophic cardiomyopathy (FHC).

# **Recommended Dilutions**

**WB** 1:500 - 1:1000

**IHC-P** 1:50 - 1:200

# **Immunogen Information**

**Gene ID**70

Swiss Prot

P68032

### **Immunogen**

A synthetic peptide corresponding to a sequence within amino acids 1-100 of human ACTC1 (NP $\_005150.1$ ).

# **Synonyms**

ACTC; ASD5; CMD1R; CMH11; LVNC4; ACTC1

# **Contact**

<b>a</b>		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
$\odot$	П	www.abclonal.com.cn


## **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

### **Storage**

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

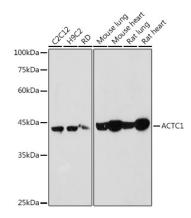


Western blot analysis of various lysates using ACTC1 Rabbit pAb (A6507) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.

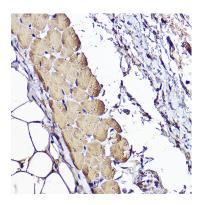

Immunohistochemistry analysis of paraffinembedded Human skeletal muscle using ACTC1 Rabbit pAb (A6507) at dilution of 1:50 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.