Leader in Biomolecular Solutions for Life Science

Phospho-CFL1-Y139 Rabbit pAb

Catalog No.: AP0330



Basic Information

Observed MW Refer to figures

Calculated MW 19kDa

Category Primary antibody

Applications IHC-P,ELISA

Cross-Reactivity Human, Mouse, Rat

Background

The protein encoded by this gene can polymerize and depolymerize F-actin and G-actin in a pH-dependent manner. Increased phosphorylation of this protein by LIM kinase aids in Rhoinduced reorganization of the actin cytoskeleton. Cofilin is a widely distributed intracellular actin-modulating protein that binds and depolymerizes filamentous F-actin and inhibits the polymerization of monomeric G-actin in a pH-dependent manner. It is involved in the translocation of actin-cofilin complex from cytoplasm to nucleus.

Recommended Dilutions

IHC-P	1:50 - 1:200
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID 1072 Swiss Prot P23528

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms

CFL; cofilin; HEL-S-15; Phospho-CFL1-Y139

Contact


6	400-999-6126
\times	cn.market@abclonal.com.cn
€	www.abclonal.com.cn

Product Information


Source Rabbit **Isotype** IgG **Purification** Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

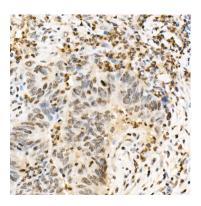

Immunohistochemistry analysis of paraffinembedded Human colon carcinoma using Phospho-CFL1-Y139 Rabbit pAb (AP0330) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

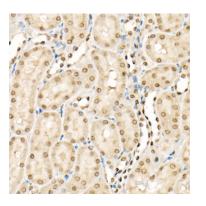

Immunohistochemistry analysis of paraffinembedded Mouse kidney using Phospho-CFL1-Y139 Rabbit pAb (AP0330) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

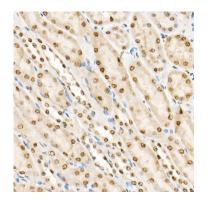

Immunohistochemistry analysis of paraffinembedded Rat kidney using Phospho-CFL1-Y139 Rabbit pAb (AP0330) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.