

ABflo® 555-conjugated Goat anti-Rabbit IgG (H+L)

Catalog No.: AS058 9 Publications

Basic Information

Observed MW

Calculated MW

Category

Secondary antibody

Applications

IF/ICC,FC

Cross-Reactivity

Wide range of species

Conjugate

ABflo® 555. Ex:557nm. Em:570nm.

Background

Secondary antibodies are affinity-purified antibodies which will work with target-specific primary antibody in the detection, sorting or purification of its specified target. Secondary antibodies offer increased versatility enabling users to use many detection systems (e.g. HRP, AP, fluorescence). They can also provide greater sensitivity through signal amplification as multiple secondary antibodies . Most commonly, secondary antibodies are generated by immunizing the host animal (different from host species of primary antibody) with a pooled population of normal immunoglobulins from the host species of primary antibody and can be further purified and modified (i.e. antibody fragmentation, label conjugation, etc.) to ensure well-characterized specificity to corresponding normal immunoglobulins.

Recommended Dilutions

IF/ICC 1:100 - 1:500

FC 1:500 - 1:2000

Immunogen Information

Gene ID Swiss Prot

Immunogen

This information is considered to be commercially sensitive.

Synonyms

Contact

<u>a</u>		400-999-6126
\bowtie		cn.market@abclonal.com.cn
\odot	T	www.abclonal.com.cn

Product Information

Source	Isotype	Purification
Goat	IgG	Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

Validation Data





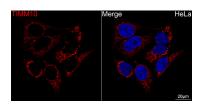

Flow cytometry: 1×10^6 K-562 cells (negative control,left) and A-431 cells (right) were surface-stained with Purified Rabbit anti-Human E-Cadherin mAb (2 μ g/mL,orange line) or secondary antibody only (blue line). Non-fluorescently stained K-562 and A-431 cells were used as blank control (red line). ABflo® 555-conjugated Goat Anti-Rabbit IgG (H+L)(AS058, 1:2000) was used as a secondary antibody.

Confocal imaging of HeLa cells using [KD Validated] TIMM10 Rabbit mAb (A24996, dilution 1:200) followed by a further incubation with ABflo® 555-conjugated Goat Anti-Rabbit IgG (H+L) (AS058, dilution 1:800). DAPI was used for nuclear staining (Blue). Objective: 100x.